MultiVIEW

SCP SCIENCE

A NEW STANDARD IN MICROWAVE DIGESTION



CORPORATE HEADQUARTERS

Phone: +1 (514) 457-0701 Fax: +1 (514) 457-4499 CANADA / USA


Phone: +1 (800) 361-6820 Fax: +1 (800) 253-5549 FRANCE

Phone: +33 (0) 1 69 18 71 17 Fax: +33 (0) 1 60 92 05 67 **GERMANY**

Phone: +49 (0) 8342-89560-61 Fax: +49 (0) 8342-89560-69 UK

Phone: +44 (800) 368 6147

WELCOME TO THE NEXT GENERATION IN MICROWAVE DIGESTION

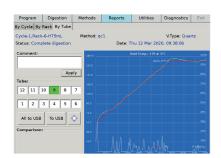



DESIGNATED INDIVIDUAL CHANNELS

- 12 individual feedback circuits control magnetron power and temperature of each sample.
- 12 magnetrons (one for each sample in a rack) provide dedicated microwave energy.
- 12 side mounted temperature detectors operate independent of volume and sample type.

USER INTERFACE

- View each sample profile on a 15" color touch screen.
- View and save temperature and microwave power (PWM) curves in real-time.
- Modify digestion parameters (hold-time, temperature-time, and ramp-time) on the fly.

PRE-ASSIGN METHODS

- Per vessel: Allow multiple methods in each rack.
- Per rack: Entire rack of the same sample type.
- On multiple racks: ideal for automation and high throughput (requires AutoLOADER).

SAFETY

- Installed parts are Teflon® coated to resist corrosive environments.
- Designed with a 13 gauge steel tunnel to protect user from exothermic reactions.
- Integrated, high speed, cooling fans exhaust fumes at 100 cfm.
- Pressure venting caps release pressure beyond 500 PSI (34 bar) to prevent vessel over pressurization.

MAXIMIZE SAMPLE THROUGHPUT

- For larger labs, increase lab productivity and upgrade to MultiVIEW AutoLOADER (sold separately).
- A patented, rotary carousel tower designed with a small footprint to save bench-top space.
- Load, stack, and process up-to 7 racks, of 1-12 samples, in a continuous cycle.

CATALOG #	DESCRIPTION	DIMENSIONS (W x H x D)	WEIGHT
010-625-001	MultiVIEW Microwave Digestion System	23" x 39" x 28" 58 x 99 x 71 cm	245 lbs 111 kg
010-625-002	MultiVIEW with AutoLOADER	48" x 40" x 28" 120 x 102 x 71 cm	395 lbs 180 kg
010-625-004	AutoLOADER, to automate the MultiVIEW and process up-to 7 racks	24" x 40" x 28" 60 x 102 x 71 cm	150 lbs 68 kg
010-625-005	MultiVIEW rack, with 75 ml Quartz vessels, calibration lines at the 10/ 25/50 ml volume	-	-
010-625-006	MultiVIEW rack, with 75 ml Teflon® vessels	-	-

CONTACT YOUR LOCAL SALES REPRESENTATIVE OR DISTRIBUTOR FOR MORE DETAILS.

CORPORATE HEADQUARTERS CANADA / USA FRANCE GERMANY UK

Phone: +1 (514) 457-0701 Fax: +1 (514) 457-4499 Phone: +1 (800) 361-6820 Fax: +1 (800) 253-5549 Phone: +33 (0) 1 69 18 71 17 Fax: +33 (0) 1 60 92 05 67 Phone: +49 (0) 8342-89560-61 Fax: +49 (0) 8342-89560-69 Phone: +44 (800) 368 6147